

AXI-PATH RECORDER ACCESSORIES

7040 AP Edentulous Outriggers

The Edentulous Outrigger System is used with either the Axi-Path II or Axi-Path III Recorder. They attach to the cross rod to add outriggers which support two elastic bands to form a sling under the patient's mandible. The anterior discluder on the Axi-Path clutch contacts the anterior edentulous ridge to hold the disposable clutch onto the edentulous ridge during the recording procedure.




7691 AP Disposable Clutches (10)

The clutches are used with the Axi-Path Recorder to attach the axis recording arms to the mandible. Each clutch comes with three different size discluders (S,M,L). The clutches are designed to break apart in four sections for easy removal after the jaw movement recording procedure is completed. They can be cut shorter by the operator so that one size fits all. The clutch is attached to the mandibular teeth with quick setting plaster.



7685 AP Axis Graph Papers (20 sets-paper)

These Graph Papers are used with either the Axi-Path II or Axi-Path III Recorder. They are printed on crack-n-peel adhesive backed paper for easy placement on the Axi-Path graph supports. The lines are marked in 1mm, 5mm and 10mm squares. The graph papers are conveniently marked with 12:00, 3:00, 6:00 and 9:00 positions for easier location of the transverse axis.

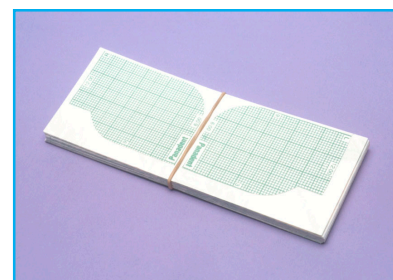

7686 AP Marking Dye (certified U.S. Government standard)

This dye is used to mark the transverse (hinge) axis on the patient's skin for axis face-bow transfer.



9531 AP Writing Leads (12 Red)

9530 AP Writing Leads (12 Black)

These are replacement leads for the writing stylus for tracing mandibular movements during the recording procedure with the Axi-Path II or Axi-Path III recorders.



7695 AP Rubber Head Bands (2)

These are replacement head bands for the Axi-Path II and Axi-Path III recorders.

