


FACE-BOW AND ACCESSORIES

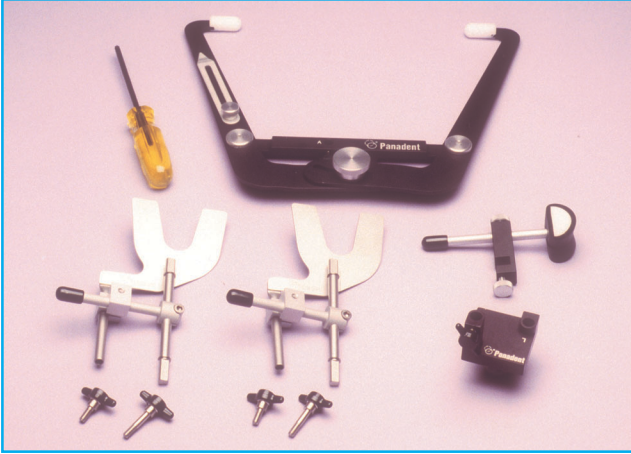

4000 FB *Pana-Mount*™ Face-Bow

The *Pana-Mount*™ Face-Bow features wrench or wrenchless detachable bite fork assemblies which allows the face-bow to be used for multiple patients without having to mount casts immediately. It has a standard 22mm nasion relator as well as an orbital pointer. Comfortable ear pieces and foam nasion pad increase patient comfort. The face-bow can be used as an average axis ear-bow as well as a true transverse axis bow (in conjunction with the Axis Mounting System and Axi-Path Recorder).

Includes:

Face-Bow, two complete Bite Fork Assemblies, Mounting Fixture, Standard Nasion Relator, T-Screws, and Hex Wrench



4080 FB Adjustable Nasion Relator

The Adjustable Nasion Relator is used in conjunction with the *Bio-Esthetic*™ *Level Gauge* to level the face-bow in the sagittal plane when the patient's head is erect. The adjustable nasion is calibrated in millimeters from nasion. There is a long and a short extension with each kit which allows the anterior reference point to be changed in a range from 25mm to 45mm (short: 25-35mm, long: 35-45mm).




4070 FB *Bio-Esthetic*™ Level Gauge

The *Bio-Esthetic*™ *Level Gauge* is used to help align the average axis face-bow horizontal so the anterior esthetic plane of the teeth can be established correctly on the articulator.


Includes:

One Magnetic Level Gauge and two self-adhesive ferrous metal discs (one for each side of face-bow)

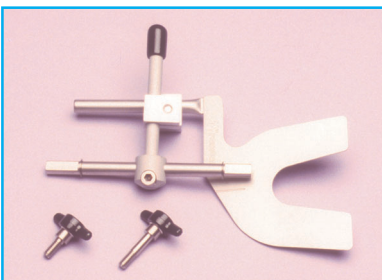

4027 FB Bite Fork Assembly (complete)

The detachable Bite Fork Assembly is used with the *Pana-Mount*™ Face-Bow and comes complete with a standard bite fork and stem assembly, as well as T-Screws. By having extra bite fork assemblies, the operator can use one face-bow for multiple patients consecutively.

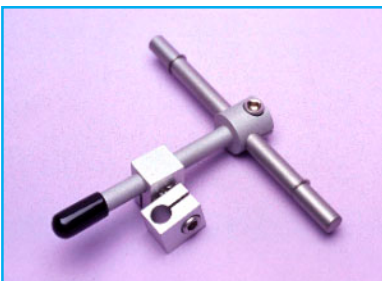

4030 FB Bite Fork Stem Assembly (only)

This fixture is for attaching the bite fork to the *Pana-Mount*™ Face-Bow. It is the same assembly as item #4027, excluding the bite fork. The T-Screws are also included.

FACE-BOW AND ACCESSORIES

4028 FB Bite Fork (regular)

The standard bite fork is used with the **Pana-Mount**™ Face-Bow and bite fork stem assembly. Constructed of all stainless steel it is designed for easy cleaning and sterilization. The mid-line mark helps for orientation in the mouth and the flat surface of the mouth piece makes for easy placement of Panadent **Bite-Tab**™ compound discs.




4037 FB Bite Fork (full denture)

The full denture bite fork is used with the **Pana-Mount**™ Face-Bow and bite fork stem assembly for transferring maxillary bite-rim and edentulous cast to the articulator. Constructed of all stainless steel it is designed for easy cleaning and sterilization. The mid-line mark helps for orientation in the mouth.



4038 FB Pedo Bite Fork (child)

The pedo bite fork is used with the **Pana-Mount**™ Face-Bow and bite fork stem assembly. Constructed of all stainless steel it is designed for easy cleaning and sterilization. The mid-line mark helps for orientation in the mouth and the flat surface of the mouth piece makes for easy placement of Panadent **Bite-Tab**™ compound discs.

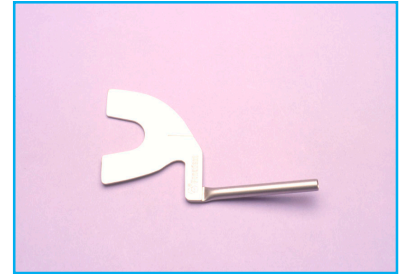

9925 FB Hex Wrench (for Face-Bow)

This wrench is used with the **Pana-Mount**™ Face-Bow for attaching the bite fork and stem assembly to the face-bow. The wrench has a large handle for easy torque and a "ball" end so the angle of insertion into the clamping toggle is not critical.



4022 FB Nasion Pad (1 set of two pieces)

The foam cushioned Nasion Pad is used with either the **Pana-Mount**™ Face-Bow or Axi-Path Recorder head frame. The two halves slide easily over the nasion relator to create a comfortable cushion for the patient's nose during face-bow or Axi-Path recording procedures.



4450 FB T-Screws for stem assembly (set of 2)

4026 FB Thumb Screw for cross bar on face bow

4025 FB T-Screw (for mounting fixture)

These T-Screws allow the Pana-Mount Face-Bow & Bite-Fork Assemblies to become a wrenchless system.

