

KOIS DENTO-FACIAL ANALYZER SYSTEM



4300 FB Kois Dento-Facial Analyzer System

The Kois Dento-Facial Analyzer System registers and transfers the patient's occlusal plane as well as tilts in the occlusal plane in three planes of space to the articulator related to an average 100mm axis-incisal distance. O-rings can be adjusted on the Vertical Indicator Rod to the patient's facial landmarks for the evaluation of facial proportions. The height of the lip commissures can be measured from the Index Tray for marking on the study casts.

Includes:

- #4305 Kois Dento-Facial Analyzer
- #4360 Kois Adjustable Platform



4305 FB Kois Dento-Facial Analyzer

The Kois Dento-Facial Analyzer is used with the Index Trays to register and transfer the patient's occlusal plane as well as tilts in the occlusal plane in three planes of space to the articulator.

Includes:

- #4322 Kois disposable Index Trays (20)




4322 FB Kois Disposable Index Trays (20)

4323 FB Kois Disposable Index Trays (50)

4324 FB Kois Disposable Index Trays (100)

The disposable index trays are used with the Kois Dento-Facial Analyzer and provides a permanent or disposable bite-fork registration record.

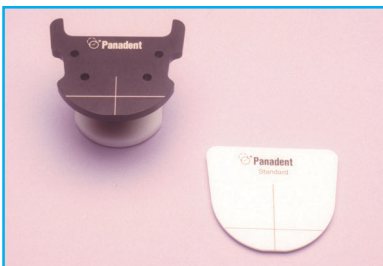

4360 FB Kois Adjustable Platform

The Kois Adjustable Platform is used with the Kois Index Trays to mount the patient's occlusal plane as well as tilts in the occlusal plane to the articulator. The platform can be adjusted in millimeters to determine incisal length. The platform is also used in conjunction with the Golden Proportion Waxing Guides.

Includes:

- # 4331 Standard Waxing Guide